


Syringe Tool, Touch-Up Paint, and Lubricant

Syringe Tool

98.0010.00

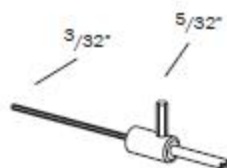

- Use to remove/replace valve cores on Uni-Flo and Tri-Flo Syringes
- Use 3/32\" to adjust syringe flow on Century II
- Use 5/32\" to remove/replace the syringe tip retaining nut

Valve Test Syringe

98.0050.01 (4)



- For making quick tests of pilot-operated valves
- Applies a static pressure of 5 to 75 psi

A-dec Silicone Lubricant

98.0090.01 0.05oz blister packs (6)

98.1120.00 Re-sealable tube 5.3oz



Touch-Up Paint

90.0351.00 Surf 4

90.0351.01 Surf 6

90.0351.02 White 2

90.0351.03 Gray 3

90.0350.09 Surf



- 1/2 fluid ounce
- Includes brush

Loctite®

060.001.00 Red #271

060.002.00 Blue #242



Service Tip

A-dec's silicone lubricant is a high quality silicone-based grease. It is ideal for lubricating A-dec O-rings. Before installing O-rings always apply a light coating of silicone lubricant. This makes the installation of the O-ring easier and prevents damage to the O-ring. A-dec silicone lubricant is available in a package of 6 as Part Number 98.0090.01.

Caution

The use of petroleum-based lubricants (jellies, creams, or lotions) causes the O-rings in the syringe to swell and/or deteriorate, resulting in syringe malfunction or failure.