


# Relay

## Time Bowl Rinse Relay

12.0913.00

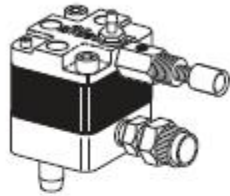

- For Cascade cuspidor
- 3/8" female Poly-Flo
- 1/8" barb
- 3/8" x 1/4" barb

## Time Delay Relay

42.0534.00




- For filter regulator manifolds and Excellence cuspidor
- 1/8" and 1/4" barb
- 9-hole gasket, 24.0137.01 (10)

## Cup Fill Relay and Cup Fill Timer Stabilizer Kit

12.0934.00 Cup fill relay

90.0456.00 Stabilizer kit

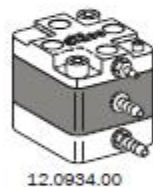

12.0934.00




90.0456.00

- For Cascade cup fill
- 1/8" and 1/4" barb

NOTE: Kit replaces the cup fill timer and improves flow adjustment.

## Relay Valve with Holdback

34.0043.00

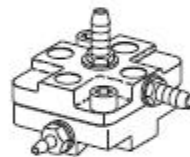

- For tooth dryer
- Non-Cascade
- 1/8" barbs (1)
- 1/4" barbs (2)

## Service Tip

Air bubbles in the water or sputtering may be caused by:

- Loose syringe handle
- Defective O-ring on the connector tube at the syringe terminal
- Leaking diaphragm in the water regulator
- Leaking diaphragm in any water relay valve