

# Shutoff and Selector

## Air-Actuated Water Shutoff

34.0031.00



- Inlet 3/8" male compression fitting
- 1/8" barb

## Air-Actuated Water Shutoff

34.0032.00 1/4" Poly-flo

34.0033.00 3/8" Poly-flo



- Inlet 3/8" male compression fitting
- 1/4" poly sleeve with 1/4" poly adapter

## Air-Actuated Water Shutoff

34.0036.00




- Inlet and outlet 3/8" female compression fittings
- For Preference dental furniture sink

## Air-Actuated Air Shutoff

34.0037.01 1/4" Poly-flo

34.0037.00 3/8" Poly-flo

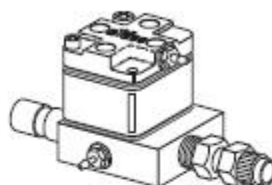

- Inlet 3/8" male compression fitting

## Manual Shutoff

026.062.00 1/2" Compression

026.158.00 1/2" Female pipe thread



026.158.00 shown

- Inlet 1/2" compression or 1/2" female pipe thread
- Outlet 3/8" compression
- Nickel plated brass


## Service Tip

Restricted water flow could be the result of a clogged screen in the manual shutoff valve. To check:

1. Turn off the master On/Off valve toggle.
2. Bleed off system pressure.
3. Turn OFF manual shutoff valve.
4. Remove water filter/regulator.
5. Check screen.

## 3-Way Selector


33.0104.05

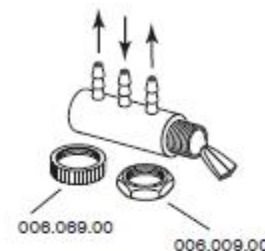

- Air applications only
- 2 position, 3-way
- 1/8" barbs
- Surf 4 plastic toggle

## Ported 3-Way


82.0305.04

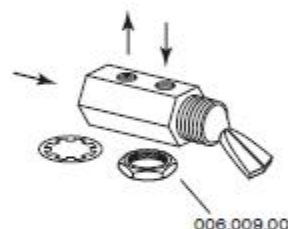

- Air or water applications
- 10 - 32 ports (3)
- Panel mount
- Surf 4 plastic toggle

## Routing/Selector 1 Function

33.0051.01

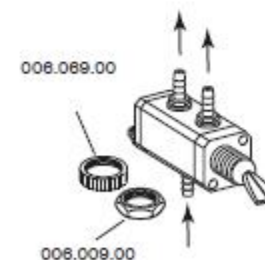

- 2 position, 2-way
- 1 function for 2 handpieces
- Drive air only
- Surf 4 plastic toggle