



Solids Collectors

Cascade Solids Collector Single HVE

For Std or 15 mm HVE

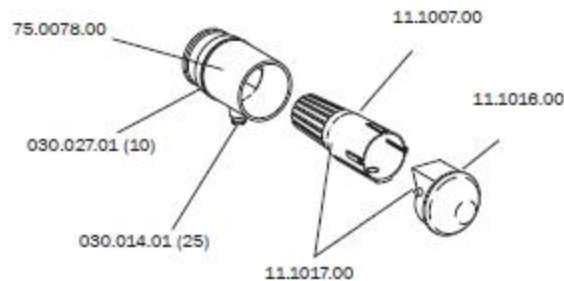

Cascade Solids Collector Dual HVE

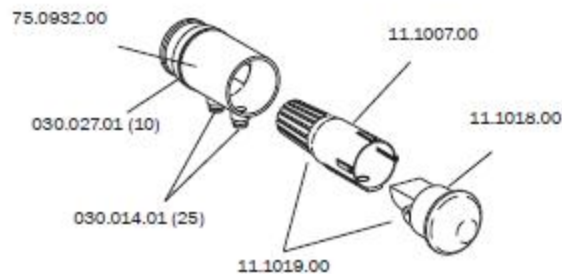

To change the A-dec dual HVE screen to a Pinnacle screen:

1. Remove the A-dec cap (11.1018.00)
Retain for later use
2. Remove the A-dec screen (11.1007.00)
3. Insert the Pinnacle screen
4. Re-insert the A-dec cap (11.1018.00)
removed in Step 1


Pinnacle Screen Kit

11.1194.00 For single HVE



- Includes (3) screens and replacement cap 11.1192.00 for single HVE
- Pinnacle screens are available from Pinnacle; order #5500-1

Cascade Solids Collector 15 mm HVE with Pinnacle Screen Kit

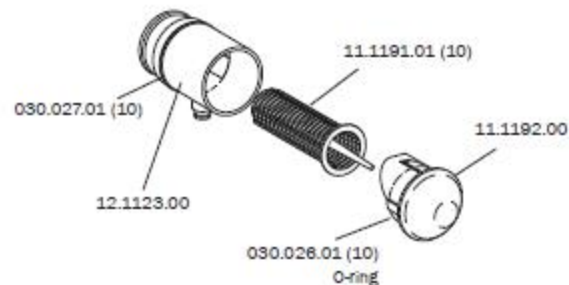

Central Vacuum Solids Collector

17.0666.00 Lid, Black

17.0670.00 Screen



- Metal canister
- 2.30" diameter screen
- Threaded

Central Vacuum Solids Collector

17.0600.02 Canister assy

028.008.00 Plug/Cap, 1/2" x 1/4"

